

## CB40-C CYTO Centrifuge



CB40-C cell smear centrifuge can be used for the liquid handling in fields of cytology, hematology and microbiology. Such as cerebrospinal fluid, pleural fluid, peritoneal fluid, synovial fluid, urine, oral cavity-cleaning, stomach cleaning, exudate, percolate, phlegm. It is mainly used for checks with small cell capacity. Ourclip rotor for cytology is compatible with popular cyto centrifuge accessories currently on the market.



### Product Advantages

- ◆Microprocessor control and DC brushless motor.
- ◆Touch panel and LCD display. RCF can be calculated automatically.
- ◆Save up to 25 user-defined programs. Featuring a separate short spin key for quick centrifugation.
- ◆Special damping structure to reduce vibration.
- ◆Centrifuge lid with soft-touch lid closure. centrifuge cannot work if the door is open and the door cannot be opened when it is working.
- ◆Several smear rotors available, it can smear various fluids to achieve the best centrifugation effect.



### Technical Parameter

Order NO	927522
Max Speed	3000r/min
Max RCF	890xg
Funnel Capacity	0.1ml-0.5ml
Rotor Capacity	12pcs
Speed Accuracy	±30r/min
Noise	≤ 60dB (A)
Power Supply	AC 220±22V 50/60Hz
Power Consumption	260W
Dimension (WxDxH)	380×500×300(mm)
Packing Dimension (WxDxH)	465×660×350 (mm)
Net Weight	19kg




### Rotor Parameters



Swing bucket cell smear rotor  
Max Capacity: 6pcs  
Max Speed: 2000r/min  
Max RCF: 545xg  
Order NO: 927489



Bucket  
Order NO: 927624

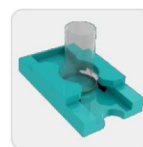

Plastic holder  
Order NO: 927657



Cell smear rotor  
Max Capacity: 6pcs  
Max Speed: 3000r/min  
Max RCF: 730xg  
Order NO: 927490



Smear clip  
Order NO: 927810




Glass slide  
Order NO: 927819



Cell smear rotor  
Max Capacity: 12pcs  
Max Speed: 3000r/min  
Max RCF: 890xg  
Order NO: 927491



Filter paper  
Order NO: 927811

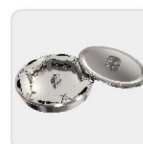

Bio-seal smear rotor  
Max Capacity: 6Plastic holder  
Max Speed: 2000r/min  
Max RCF: 405xg  
Order NO: 927455