



CCL Series Climate Chamber with Light / Plant Growth Chamber

With light: 10~50°C | Without light: 0~50°C | 150L - 1500L

LAB1ST CCL Series Climate Chamber provides a precisely controlled environment tailored for advanced plant cultivation, seed germination, and photobiological research. It integrates high-precision temperature control with energy-efficient LED lighting to simulate ideal growth conditions.

Engineered for demanding applications in agricultural research, forestry, and biotechnology, the CCL series ensures consistent, reproducible results for your most critical experiments.



Features


> User-Centered & Intelligent Operation

- **Intuitive Touchscreen Interface:** Features a large LCD screen that displays multiple parameters simultaneously with a clear, menu-driven interface for straightforward setup, monitoring, and operation.
- **30-Stage Programmable Control:** Simplifies complex experimental protocols. Each program supports flexible timing from **1** minute to **99** hours and **59** minutes, enabling true walk-away automation.
- **Versatile Simulation Capabilities:** Can simulate natural day/night temperature variations and multi-directional light sources.

> Advanced, Energy-Efficient LED Lighting System

- **Precision Spectrum & Cool Operation:** Utilizes energy-saving LED lighting that provides stable, full-spectrum illumination ideal for photosynthesis without generating excess heat, ensuring consistent chamber temperature and optimal plant growth.
- **Superior Performance:** Compared to traditional lighting, our LED system offers precise wavelength control, an extended lifespan exceeding **30,000** hours, and reduces energy consumption by approximately **80%**.
- **Rapid Response & Uniformity:** Optional planar lighting configuration ensures even light distribution for consistent absorption by specimens, maintaining experimental reproducibility.

>LED Spectrum Diagram

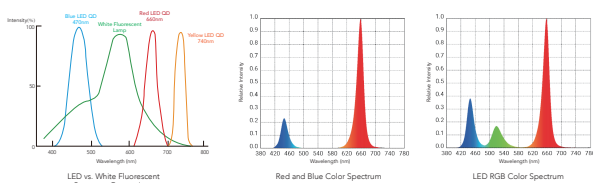

> Hygienic & Flexible Chamber Design

- **Easy-Clean Construction:** Features a mirror-finish stainless steel interior with smooth rounded corners. Internal shelves or partitions can be removed without tools, facilitating thorough chamber disinfection and cleaning.
- **Gentle Airflow Design:** The optimized air circulation system provides uniform conditions without excessive wind speed that could disturb or damage delicate plant seedlings.

> Unmatched Reliability for Continuous Research

- **Dual Redundant Compressor System:** A key differentiating feature. Two imported compressors operate in an automatic, alternating cycle to ensure uninterrupted, long-term operation—overcoming a common limitation in standard growth chambers.
- **Power Failure Protection:** User-defined parameters are automatically saved in non-volatile memory, ensuring programs resume automatically after a power outage.

> Comprehensive Safety & Protection


- **Independent Over-Temperature Alarm:** Automatically cuts off power and triggers alerts if temperature limits are exceeded, safeguarding both your experiments and the equipment.
- **Built-In Diagnostics:** The system performs self-checks and displays clear fault messages on the screen for immediate identification and troubleshooting of operational issues.
- **Audible & Visual Alarms:** Includes alerts for high temperature and over-temperature conditions.

> Eco-Friendly Design

- **Fluorine-Free Construction:** Manufactured with environmentally friendly materials, supporting a greener and healthier approach to research.



Model	CCL-150	CCL-250	CCL-300	CCL-450	CCL-800	CCL-1000	CCL-1500
Temperature Control Range	With light: 10~50°C Without light: 0~50°C						
Temperature Resolution	0.1°C						
Temperature Fluctuation	±1°C						
Light Intensity Range	0 ~ 18,000 Lux	0 ~ 18,000 Lux	0 ~ 28,000 Lux	0 ~ 35,000 Lux	0 ~ 42,000 Lux	0 ~ 50,000 Lux	0 ~ 55,000 Lux
Light Intensity	6-level adjustable						
Light Source / Illumination Method	LED, door-mounted lighting		LED / Three-side lighting		Fluorescent lamp, partition-type lighting		
Power Supply	AC220V 50HZ						
Input Power (W)	1500	1750	1700	2000	2700	2800	4300
Chamber Volume (L)	150	250	300	450	800	1000	1500
Interior Dimensions (WxDxH mm)	550×400×670	600×610×830	520×550×1140	700×550×1140	800×590×1650	1050×590×1650	1550×590×1650
Exterior Dimensions (WxDxH mm)	690×800×1410	760×815×1550	830×850×1850	950×850×1850	1475×890×1780	1410×890×1950	1570×1475×2050
Program Control Function	Independent setting of temperature, and light intensity. Capable of setting 30 program segments, with each segment's time range adjustable from 1~99 hours 59 minutes.						
Operating Environment Temperature	+5~30°C						
Continuous Operation Time	Capable of long-term continuous operation (two sets of international fully sealed compressors automatically switch in rotation).						

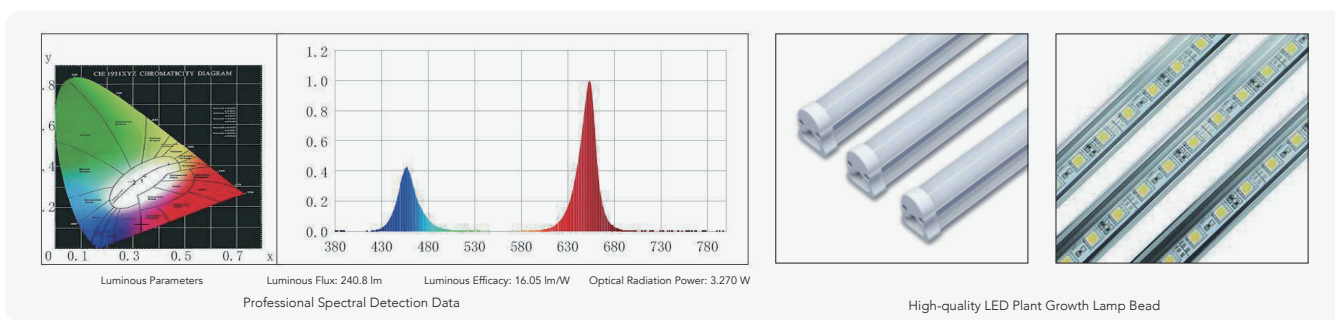

Options

- RS485/232 Interface and Communication Software**
- Embedded Printer**
- USB Data Storage**
- CO₂ Inlet Port** (for promoting plant growth)
- CO₂ Controller** (with Infrared CO₂ Sensor)
- Illumination Control System** (Stepless Dimming)

Note:

CO₂ Control System (Optional)

- For advanced plant physiology studies, an optional infrared CO₂ sensor and control system is available. It provides rapid, accurate response to concentration changes (range: 0-5000ppm or 0-20%) unaffected by temperature or humidity, ensuring precise atmospheric control.

Labfirst Scientific Instruments (Shanghai) Co., Ltd.
No.248 Guanghua Road, Minhang District, Shanghai, China

CALPHA INDUSTRIES INC.
10 Corporate Park Ste 330, Irvine, California 92606, USA

LSI PROCESS EQUIPMENT CO., LTD.
701 West Georgia Street, Suite 1500, Vancouver, British Columbia, V7Y 1C6, Canada

Tel: +86-13524020331 / +1-844-452-2178

sales@lab1st.com
www.lab1st.com
eshop.lab1st.com